

APPLICATION NOTE



Cannabinoid Purification and Scale-Up on Hipersep[®] Pilot

LOGAN CLINE (1), OLIVIA COOLEY (1), ALPESH PATEL (2) AND JIN SEOK HUR (2),
(1) Extract Labs™ 3620 WALNUT STREET, BOULDER, CO 80301 USA
(2) NOVASEP, LLC., 23 CREEK CIRCLE, BOOTHWYN, PA 19061 USA





Figure 1. Varicol® (advanced SMB) and Hipersep® with 99.99% solvent recycling for MTon Scale purification of natural extract for pharmaceutical application.

Introduction

For the outstanding efficiency and scalability, high performance liquid chromatography (HPLC) has been widely used for purification of high quality APIs and intermediates in pharmaceutical industry. The application has been broadened to purification of natural extracts, which are often very complex mixtures.

This application note presents the preparative HPLC purification of cannabinoids. The process was developed using a bench-top HPLC instrument and then directly scaled-up to Hipersep® Pilot with Prochrom® LC150 (150mm i.d.) for manufacturing.

Acknowledgement

We thank Extract Labs™ for sharing the data and Nouryon for providing the bulk C18.



Purification of Cannabinoids

Cannabinoid oil extracted from cannabis plants contains various cannabinoids such as Δ^9 -tetrahydrocannabinol (Δ^9 -THC), which is a well-known psychoactive compound, and cannabidiol (CBD). Although one of the main objectives is to remove THC and pesticides, the demand of purifying rare cannabinoids such as cannabigerol (CBG) cannabinol (CBN) and cannabichromene (CBC) has been considerably increased.

Since the cannabinoids have very similar structures and similar chemical and physical properties, purifying them with high purity and yield is certainly a challenging task. In addition, robust scale-up of the purification process from few gram to Mton scale requires another level of expertise not only on the purification process but also the engineering requirement.

In order to achieve the target products with a high purity at a preparative scale, cannabinoid oil was first obtained from the plant by extraction followed by several unit operations. Then, the purification condition was developed and optimized at bench scale using an analytical size column packed with C18 stationary phase. The purification condition is described below.

Stationary phase:	Kromasil C18, 13 μ m, 100Å
Column size:	4.6 i.d. x 250mm
Mobile phase:	Mixture of Methanol and Water
Flow rate:	0.5 to 1mL/min
Temperature:	Room temperature to 25°C
Detection:	UV at 220 to 250nm

After the peak identification and optimization, the purification process was scaled up to Prochrom® LC150 using the same C18, which is available in bulk.

Process Scale-Up on Hipersep[®] Pilot

For predictable scale-up from the bench-top instrument to the production instrument, both the column diameter and length were increased to 150mm i.d. and 350mm length in order to maximize the production rate. Then, the flow rate was adjusted to maintain the run time per injection considering the changes on the column diameter and length. Overall, the column volume was increase from 4.1mL to 6.2L.

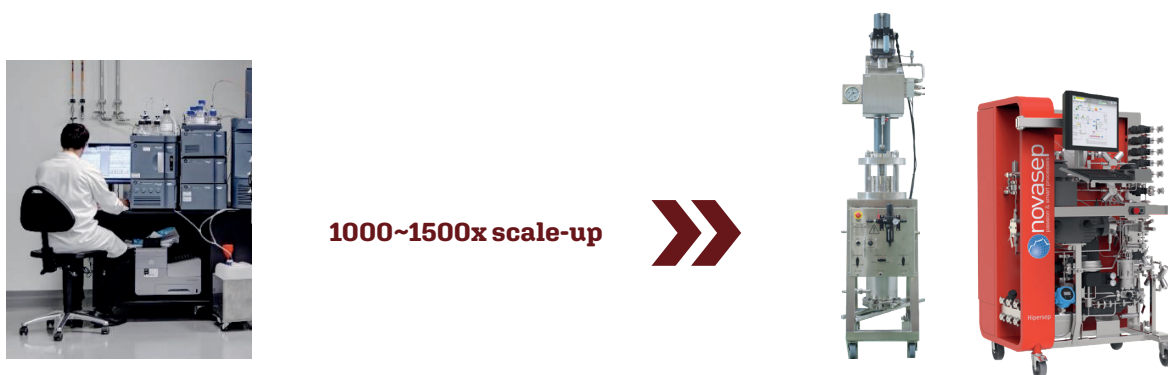


Figure 2. Approximately 1000~1500 times scale-up from analytical size column to Prochrom[®] LC150 with Hipersep[®] Pilot.

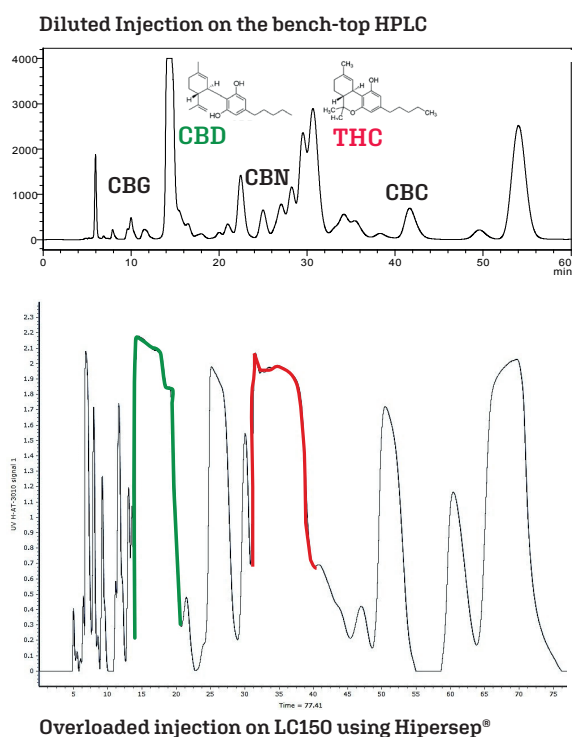


Figure 3. Comparison of the preparative chromatogram from Hipersep[®] Pilot and Prochrom[®] LC150 with diluted injection on the bench-top.

As shown in Figure 3, the purification was linearly scalable so the elution of each cannabinoid could be easily identified on the preparative chromatograms from Hipersep[®] Pilot with Prochrom[®] LC150 compared to the analytical chromatogram from the bench-top HPLC using the 4.6 mm i.d. column. In addition to CBD and THC, the process allows one to purify several rare cannabinoids with a high purity as well.

Conclusion

The purification process was successfully scaled-up from an analytical size column to Prochrom[®] LC150 using Hipersep[®] Pilot. Since the process is confirmed, it can be further linearly scaled up to larger size columns such as Prochrom[®] LC1000. It should be noted that Novasep has a patent pending for another CBD purification technique.

Prochrom Hipersep®


Skids and columns

	HIPERSEP®							
	Pilot	XS	S	M	L	XL	XXL	Above
Productivity range (kg/day) in Batch mode	0,2-4	2-7	2-7	7-33	17-57	29-106	53-165	Customization
Productivity range (kg/day) in Cyclojet	8-12	14-21	14-21	66-99	114-171	212-318	330-495	
Typical Prochrom® column diameter (mm)	50/80/110/150	150/200	150/200	300/450	450/590	590/800	800/1000	
Skid Flowrate (L/h)	6*-90	60-200	90-200	150-500	300-1000	500-1600	1000-2500	
Gradient accuracy (%)	< 1							
Max pressure (bar)	100						70	
Explosion proof	Atex Zone II Class 1 Division 2							
Certifications / Equipment standards	CODAP/ASME for columns CE/UL for systems							
Configurable options	High pressure injection module, Eluent & feed filters, Heat exchanger, Backflush configuration valves, pH conductivity, Nitrogen flush for product recovery, Up to 10 collections, Solvent recycling, Cyclojet process							


*Gradient accuracy <2% below 10L/h




Hipersep® Pilot

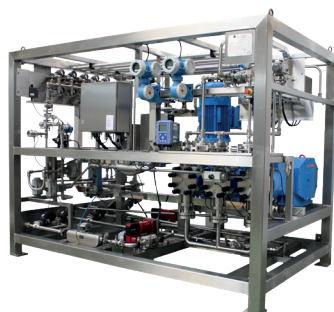

Hipersep® XXL

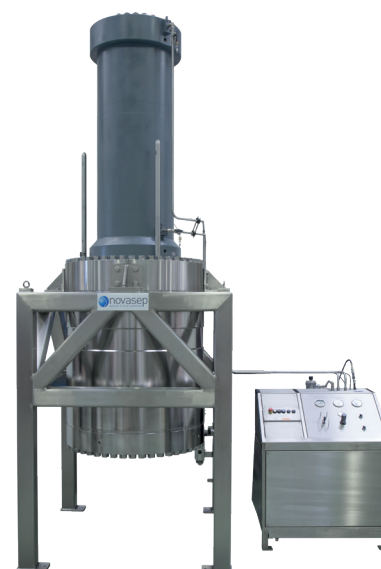


Prochrom® Hipersep® LC1000

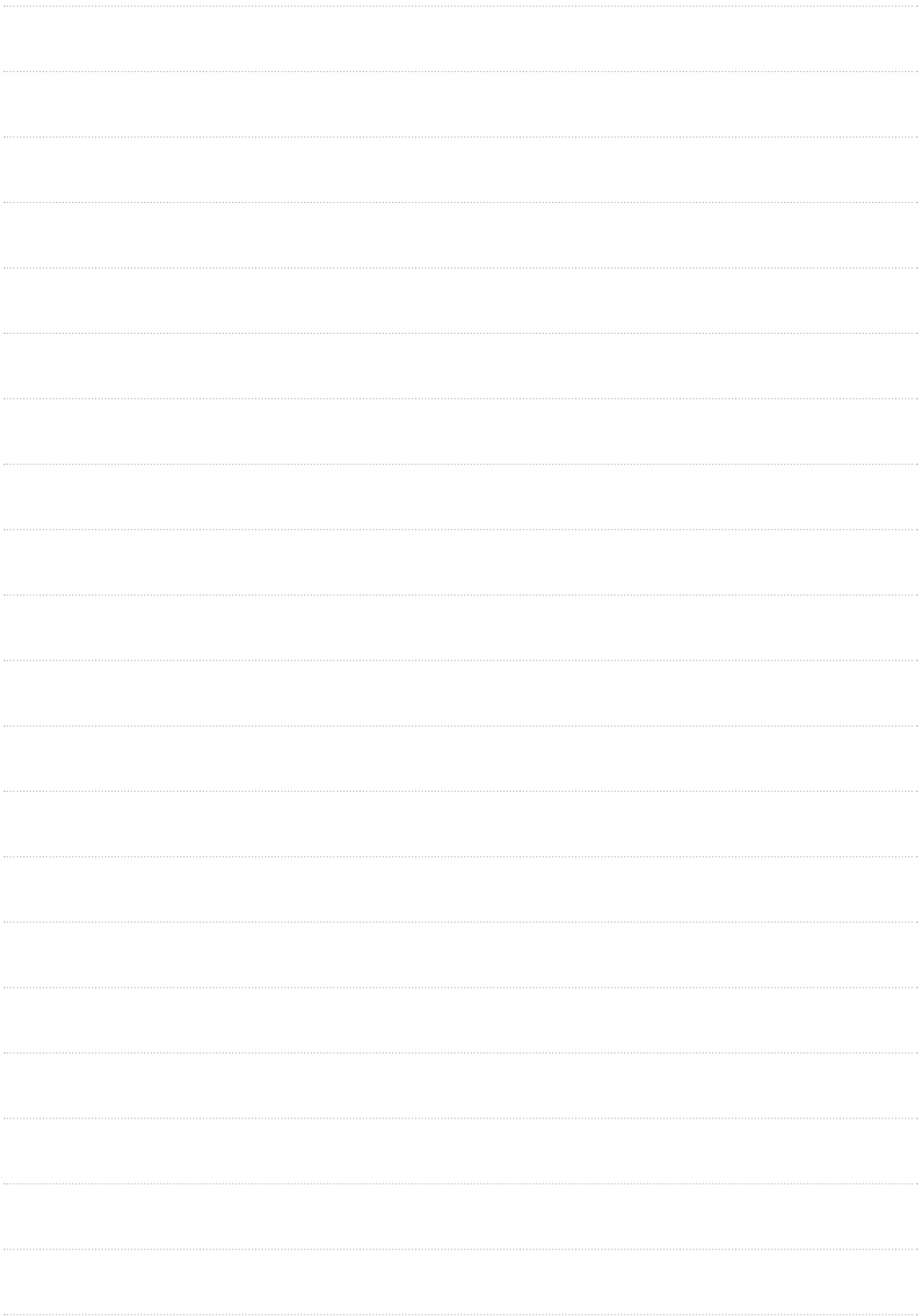

About Novasep

Novasep is a leading provider of services in the field of molecule production and purification for the life science and chemical industries, based on an unrivalled pool of specialized technologies.

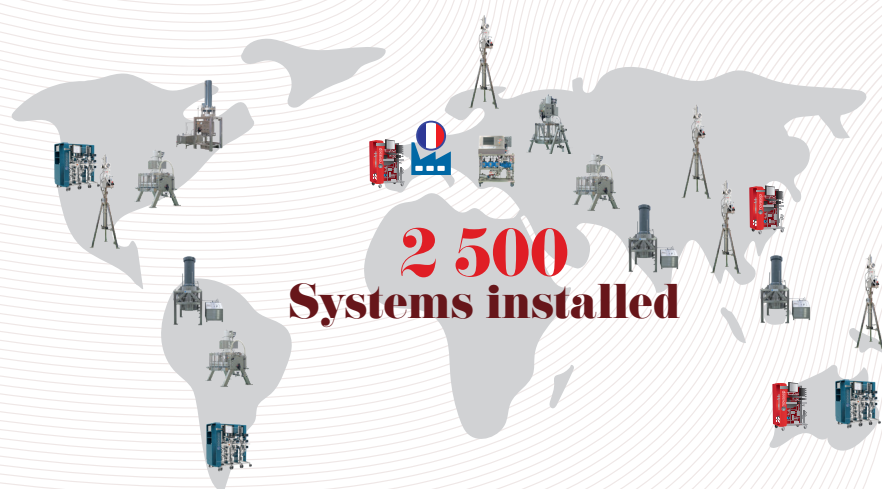

Services and technologies
for the life science industries

www.novasep.com

Contact our experts:
novasep@novasep.com



Follow us on  

EUROPE
Novasep Process
Site Eiffel
81 boulevard de la Moselle - BP 50
54340 Pompey - FRANCE
+33 383 497 000

AMERICAS
Novasep LLC
23 Creek Circle
Boothwyn, PA
19061 - USA
+1 484 361 6000

ASIA
Novasep Asia
Building No.4,
Lane 1690, Zhangheng Road,
Zhangjiang Hi-Tech Park,
Shanghai 201203 - CHINA
+86 21 6045 1600

INDIA
Novasep India
Four Square Elegant, Block A,
Flat 110,
Madhavapuri Hills, Road 5,
Hyderabad- 500050 - INDIA
+91 22 41515151